Inspected Dimensions:
Folded width: 5-1/2" (14 cm)
Folded length: 7-1/2" (19 cm)

ARROW®

Peripherally Inserted Trimmable Central Venous Catheter Product with Placement Wire

Safety and Efficacy Considerations:

Do not use if package has been previously opened or damaged. **Warning: Prior to use read all package insert warnings, precautions, and instructions. Failure to do so may result in severe patient injury or death.**

Do not alter the catheter or any other kit/set component during insertion, use or removal except as instructed.

Procedure must be performed by trained personnel well versed in anatomical landmarks, safe technique and potential complications.

Warning: Do not place the catheter into or allow it to remain in the right atrium or right ventricle (refer to Fig. 1)

practice^{1,2,4,5,7,17,18} that may lead to cardiac perforation and tamponade.^{1,2,4,5,17,18} Although cardiac tamponade secondary to pericardial effusion is uncommon, there is a high mortality rate associated with it.²¹ Practitioners placing central venous catheters must be aware of this potentially fatal complication before advancing the catheter too far relative to patient size. No particular route or catheter type is exempt from this potentially fatal complication.¹⁸ The actual position of the tip of the indwelling catheter should be confirmed by x-ray after insertion.^{1,2,5,17,18,22} Central venous catheters should be placed in the superior vena cava^{1,2,4,5,7,17,23} above its junction with the right atrium and parallel to the vessel wall^{10,23} and its distal tip positioned at a level above either the azygos vein or the carina of the trachea, whichever is better visualized. Central venous catheters should not be placed in the right atrium unless specifically required for special relatively short term procedures such as aspiration of air emboli during neurosurgery. Such procedures are, nevertheless, risk prone and should be closely monitored and controlled.

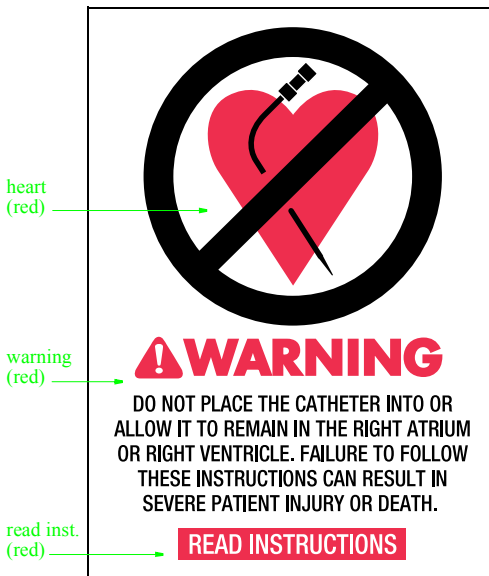
Fig. 1

Cardiac Tamponade: It has been documented by many authors that placement of indwelling catheters in the right atrium is a dangerous

Indications for Use:

A Peripherally Inserted Central Venous Catheter permits venous access to the central circulation through a peripheral vein. It offers an alternative method of intravenous therapy for select adult and pediatric patients.

Contraindications:

None known.

Warnings and Precautions:*

- Warning: Sterile, Single use: Do not reuse, reprocess or resterilize. Reuse of device creates a potential risk of serious injury and/or infection which may lead to death.**
- Warning: Do not place the catheter into or allow it to remain in the right atrium or right ventricle. Central vein catheters should be positioned so that the distal tip of the catheter is in the superior vena cava (SVC) above the junction of the SVC and the right atrium and lies parallel to the vessel wall.**
- Warning: Practitioners must be aware of complications associated with central vein**

catheters including cardiac tamponade secondary to vessel wall, atrial or ventricular perforation, pleural and mediastinal injuries, air embolism, catheter embolism, thoracic duct laceration, bacteremia, septicemia, thrombosis, inadvertent arterial puncture, nerve damage, hematoma, hemorrhage and dysrhythmias.

4. **Warning:** Do not apply excessive force in placing or removing catheter. If placement or withdrawal cannot be easily accomplished, an x-ray should be obtained and further consultation requested.
5. **Warning:** The practitioner must be aware of potential air embolism associated with leaving open needles or catheters in central venous puncture sites or as a consequence of inadvertent disconnects. To lessen the risk of air embolism, only securely tightened Luer-Lock connections should be used with this device. Follow hospital/agency protocol to guard against air embolism for all catheter maintenance.
6. **Warning:** The practitioner must be aware of clinical conditions that may limit the use of PIC catheters such as: dermatitis, cellulitis and burns at or about the insertion site, previous ipsilateral venous thrombosis, radiation therapy, contractures, mastectomy, and potential use for AV fistula.
7. **Warning:** To minimize the risk of damage to the catheter if the placement wire is removed prior to catheter insertion, do not attempt to reinsert the placement wire.
8. **Warning:** Due to the risk of exposure to HIV (Human Immunodeficiency Virus) or other blood borne pathogens, health care workers should routinely use universal blood and body-fluid precautions in the care of all patients.
9. **Precaution:** Indwelling catheters should be routinely inspected for desired flow rate, security of dressing, correct catheter position and for secure Luer-Lock connection. Use centimeter markings to identify if the catheter position has changed.
10. **Precaution:** Only x-ray examination of the catheter placement can ensure that the catheter tip has not entered the heart or no longer lies parallel to the vessel wall. If catheter position has changed, immediately perform x-ray examination to confirm catheter tip position.
11. **Precaution:** Alcohol and acetone can weaken the structure of polyurethane materials. Check ingredients of prep sprays and swabs for acetone and alcohol content.

Acetone: Do not use acetone on catheter surface.

Acetone may be applied to skin but must be allowed to dry completely prior to applying dressing.

Alcohol: Do not use alcohol to soak catheter surface or to restore catheter patency. Care should be taken when instilling drugs containing high concentration of alcohol. Always allow alcohol to dry completely prior to applying dressing.

These agents may also weaken the adhesive bond between the StatLock® Anchoring Device and skin.

12. **Precaution:** Some disinfectants used at the catheter insertion site contain solvents, which can attack the catheter material. Assure insertion site is dry before dressing.
13. **Precaution:** The indications for use in children are the same as adults; however, insertion techniques are often modified according to the age and size of a child. If the practitioner is inexperienced in utilizing this product in a child, appropriate consultation should be sought.
14. **Precaution:** When catheter is placed with distal end located in a vessel proximal to the SVC (midclavicular or midline), practitioners must be aware of the risk of chemically induced thrombophlebitis.¹¹ When administering > 10% glucose solution/total parenteral nutrition, any medication with an osmolarity above 600 mOsm/liter, medications with an osmolarity between 400 and 600 mOsm/liter for more than 4 weeks, or any medication known to be irritating to vessels proximal to the vena cava, the catheter tip must be located in the central circulation.¹⁹
15. **Precaution:** To minimize the risk of pressure induced damage to catheter, do not expose to pressures above 50 psi. Common sources of potentially high pressure include: syringes smaller than 10 cc used to irrigate or declot an occluded catheter (a fluid filled 1 cc syringe can exceed 300 psi⁶), certain radiographic procedures, and infusion pumps with occlusion pressure limits above 50 psi.
16. **Precaution:** The StatLock® Anchoring Device should be routinely inspected for adherence to the skin and connection to the catheter. At the same time, the catheter position should be checked using the centimeter marks on the catheter body.

A Suggested Procedure: Use sterile technique.

1. Verify physician order. Physician order must include post-placement chest x-ray.
2. Explain procedure to patient. Have informed consent signed as required.

- Measure and record upper arm circumference.
- Utilizing tourniquet as necessary, identify appropriate vein for insertion.
Adult: PIC catheters are typically inserted into basilic, median cubital or cephalic veins (refer to Fig. 2).

Pediatric: PIC catheters are typically inserted into basilic, cephalic, scalp veins or saphenous veins.⁹ The most suitable vein depends on size and maturity of child.

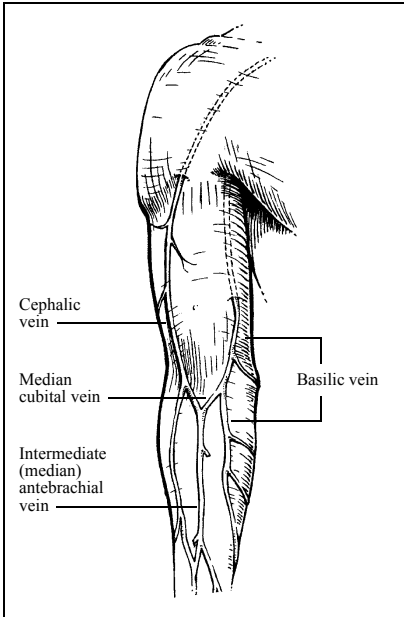
Fig. 2

- Release tourniquet leaving in place beneath arm.
- Measure patient to determine length of catheter required to place distal tip in SVC.
Adult: Extend arm approximately 45 to 90 degrees from trunk. Measure distance from insertion site along presumed anatomical course of vessel to be catheterized. Catheter tip should lie in distal one-third of SVC above right atrium and parallel to SVC wall. If a StatLock® Anchoring Device will be used, add 1 to 1-1/2 inches (2.5 to 3.8 cm) to catheter measurement.
Pediatric: Measure distance from insertion site along presumed anatomical course of vessel to be catheterized. Catheter tip should lie in distal one-third of SVC above right atrium and parallel to

SVC wall. When insertion site is saphenous vein, catheter tip should lie in distal one-third of inferior vena cava (IVC) below right atrium and parallel to IVC wall. If a StatLock® Anchoring Device will be used, add 1 to 1-1/2 inches (2.5 to 3.8 cm) to catheter measurement.

- Position patient for insertion.
Adult: Position patient as appropriate for insertion site. Extend arm laterally 45 to 90 degrees from trunk.
Pediatric: Position patient as appropriate for insertion site.
- Perform surgical scrub. Dress in protective clothing (mask, goggles, sterile gown, sterile gloves, hair cover).
- Prepare venipuncture site utilizing appropriate antiseptics.
- Perform a skin wheal with a local anesthetic as needed. In kits where provided, a SharpsAway® disposal cup is used for disposal of needles. Push needles into foam after use. Discard entire cup at completion of procedure. **Precaution: Do not reuse needles after they have been placed into disposal cup. Particulate matter may adhere to needle tip.**
- Prepare all equipment. Remove catheter tip protector. Retract placement wire approximately 1 to 1-1/2 inches (2.5 - 3.8 cm) from where catheter is to be cut. **Warning: Do not cut placement wire.** Kink proximal end of placement wire at connector with sideport. This minimizes risk of placement wire extending beyond distal tip of catheter during insertion.^{14,15} Peel back contamination guard exposing catheter portion to be trimmed. Using trimming device, cut catheter straight across (90° to catheter cross-section) to maintain a blunt tip.
- Attach syringe to sidearm and flush distal lumen with sterile saline solution. Leave syringe in place. Place catheter on sterile field.
- Reapply tourniquet and replace sterile gloves.
- Drape site for maximal barrier protection.
- Using peel-away sheath over needle, perform venipuncture. When flashback is established, advance needle and peel-away sheath as a unit until sheath is adequately within vessel.
- Release tourniquet.
- Hold sheath steady and remove needle. **Warning: To minimize the risk of possible sheath embolism, do not reinsert needle into sheath.** Check for pulsatile flow as an indicator of inadvertent arterial puncture. **Precaution: The color of blood is not always a reliable indicator of venous access.¹³**

18. Grasp distal end of contamination guard, pull back to expose catheter tip, and advance catheter slowly through peel-away sheath. The contamination guard will peel off catheter as they are advanced simultaneously.
19. When catheter tip reaches level of axillary line, position patient to decrease possibility of catheter tip malposition into internal jugular vein.
Adult: Ask patient to turn their head toward catheterized arm and drop their chin to shoulder.²⁰
Pediatric: Instruct patient as appropriate. If resistance is met while advancing catheter, retract and/or gently flush while advancing.
20. Stop advancing catheter 2 inches (5 cm) before reaching pre-established insertion length.
21. Withdraw peel-away sheath until free from venipuncture site.
22. Grasp tabs of peel-away sheath and pull apart, away from catheter, until sheath splits down entire length.
23. Advance catheter to final indwelling position.
24. Remove Luer-Lock sidearm assembly and placement wire as a unit
25. Check catheter placement with syringe by aspirating through distal lumen until free flow of venous blood is observed. **Precaution: The color of blood is not always a reliable indicator of venous access.**¹³
26. Flush lumen with sufficient volume of solution to completely clear blood.
27. Connect pigtail to appropriate Luer-Lock line as required. Alternately, port may be “locked” through injection cap using standard hospital/agency protocol. **Warning: This product contains slide clamp which may be inadvertently removed and potentially aspirated, by children or confused adults. In such situations, practitioner should remove clamp when not in use.** Slide clamp is provided on pigtail to occlude flow through lumen during line and injection cap changes. **Precaution: To minimize the risk of damage to pigtail from excessive pressure, clamp must be opened prior to infusing through lumen.**

Arrow UserGard® Needle-Free Injection Hub (where provided)

Instructions for Use:

- Attach Luer end of UserGard® hub to syringe.
- Prepare injection cap with alcohol or povidone iodine per standard hospital/agency protocol.
- Remove red dust cap from hub.
- Press UserGard® hub onto injection cap and twist to lock on pin (refer to Fig. 3).

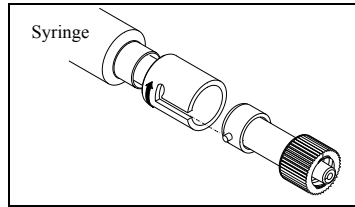
Fig. 3

- Inject or withdraw fluid as required.
 - Disengage UserGard® hub from injection cap and discard. **Warning: To minimize the risk of possible air embolism, do not leave UserGard® hub connected to injection cap.** Single use only.
28. Cleanse insertion site per hospital/agency protocol.
 29. Secure catheter. Where provided, a catheter clamp, fastener, StatLock® Anchoring Device or Steri-Strips®† may be used.

StatLock® Anchoring Device Instructions:

Initial Application:

- Cleanse and prep anticipated dressing site per hospital/agency protocol. Skin prep should be applied to coat skin and maximize StatLock® adherence. Allow to dry thoroughly. The anchor pad will be placed so center of pad is within 1 to 1-1/2 inches (2.5 to 3.8 cm) of catheter insertion site.
- The catheter can be secured to StatLock® by using the primary suture hub. **Precaution: Minimize catheter manipulation throughout this procedure to maintain proper catheter tip position. Provide strain relief placement of catheter under dressing.**

Primary Suture Hub Securement:

From earlier measurements, after catheter is properly positioned the primary suture hub is within 1 to 1-1/2 inches (2.5 to 3.8 cm) from insertion site. Place suture hub wings over StatLock® posts and press down (refer to Fig. 4). Snap StatLock® retainer wings to closed position to secure suture hub (refer to Fig. 5).

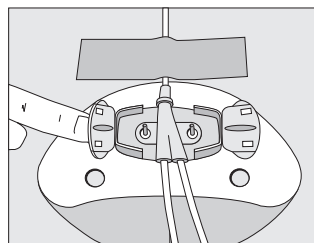
Fig. 4

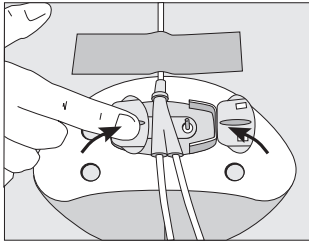
Fig. 5

- Remove paper backing from one half of StatLock® anchoring pad and press onto dry, prepared skin. Repeat process for other half of StatLock®.
- Complete sterile insertion site dressing according to established hospital/agency protocol.
- Document StatLock®/dressing application on patient's chart.
- Replace StatLock®/dressing per hospital/agency protocol. StatLock® anchoring device should be replaced at least every 7 days to ensure maximum adherence.

Alternate Technique:

Apply Steri-Strips®† and sterile dressing or apply suture and sterile dressing according to hospital/agency protocol. **Precaution: Do not suture directly to outside diameter of catheter to minimize the risk of cutting or damaging the catheter or impeding catheter flow.**

StatLock® Anchoring Device Removal:

- Remove dressing and sutures (if applicable).
- Open StatLock® retainer wings and remove catheter from StatLock® posts.
- Apply alcohol swab to StatLock® adhesive and gently lift pad off of skin.
- Document StatLock® removal procedure on patient's chart.

30. Apply PICC label to dressing.
31. Obtain chest x-ray immediately to verify tip placement. **Precaution: X-ray exam must show the catheter located in the right side of the mediastinum in the SVC above its junction with the right atrium and parallel to the vessel wall^{10,23} and its distal tip positioned at a level above either the azygos vein or the carina of the trachea, whichever is better visualized.** If catheter tip is malpositioned, reposition and re-verify. **Precaution: If difficulty is encountered in visualizing the catheter tip with x-ray, a small amount of radiopaque dye may be injected to assist in locating the tip.¹⁴**

32. Complete adhesive chart label and affix to patient's chart. Document assessments and insertion procedure.
33. Check catheter tip position routinely using centimeter marks on catheter.

Heparinization:

1. A variety of "locking" solution concentrations may be utilized to maintain patency of catheter. The amount of heparin used, if any, and frequency of flushing depends on physician preference, hospital/agency protocol, and patient condition.^{3,8}
2. The volume of heparin solution should be equal to, or slightly more than, the volume of the lumen that is being locked. Catheter priming volume is printed on product packaging.

Catheter Removal Procedure:

1. Remove dressing and sutures (if applicable). **Precaution: To minimize the risk of cutting catheter, do not use scissors to remove dressing.**
2. Open StatLock® retainer wings and remove catheter from StatLock® posts (if applicable).
3. Remove catheter by slowly pulling it parallel to skin. **Precaution: To minimize the risk of catheter breakage, do not exert excessive force if difficulty is encountered upon removal.** If resistance is met, apply heat for 20-30 minutes to area.^{16,25} Gently begin pulling catheter parallel to skin. If further difficulty is encountered, obtain an x-ray and consult physician.
4. Upon removal of catheter, measure and inspect to ensure that entire catheter length has been removed.
5. Apply alcohol swab to StatLock® adhesive and gently lift pad off of skin (if applicable).
6. Dress insertion site.
7. Document catheter removal and StatLock® removal (if applicable) removal procedures on patient's chart.

References:

1. Bar-Joseph G, Galvis AG. Perforation of the heart by central venous catheters in infants: guidelines to diagnosis and management. *J Pediatr Surg.* 1983;18:284-287.
2. Blitt CD, ed. *Monitoring in Anesthesia and Critical Care Medicine.* Central venous pressure monitoring. New York, NY: Churchill Livingstone; 1985:121-165.
3. Bonn J. Venous access: Peripherally inserted central catheters. 1994 Program, SCIVR '94, 19th Annual Scientific Meeting. San Diego, CA
4. Brandt RL, Foley WJ, Fink GH, Regan WJ. Mechanism of perforation of the heart with production of hydropericardium by a venous

- catheter and its prevention. *Am J Surg.* 1970;119:311-316.
5. Collier PE, Ryan JJ, Diamond DL. Cardiac tamponade from central venous catheters – report of a case and review of the English literature. *Angiology.* 1984;35:595-600.
 6. Conn C. The importance of syringe size when using implanted vascular access devices. *J Vasc Access Nurs.* 1993;3:11-18.
 7. Curelaru J, Linder LE, Gustavsson B. Displacement of catheters inserted through internal jugular veins with neck flexion and extension. *Intens Care Med.* 1980;6:179-183.
 8. Danek GD, Noris EM. Pediatric IV catheters: Efficacy of saline flush. *Pediatric Nursing.* 1992;18:111-113.
 9. Frey MA. Pediatric peripherally inserted central catheter program report – a summary of 4,536 catheter days. *JIN.* 1995;18:280-291.
 10. Iberti TJ, Katz LB, Reiner MA, Brownie T, Kwon KB. Hydrothorax as a late complication of central venous indwelling catheters. *Surgery.* 1983;842-846.
 11. Intravenous Nurses Society (INS). Midline and midclavicular catheters. *JIN.* July/August 1997;20:175-178.
 12. James L, Bledsoe L, Hadaway LC. A retrospective look at tip location and complications of peripherally inserted central catheter lines. *JIN.* 1993;16:104-109.
 13. Jobs DR, Schwartz AJ, Greenhow DE, Stephenson LW, Ellison N. Safer jugular vein cannulation: recognition of arterial punctures and preferential use of the external jugular route. *Anesthesiology.* 1983;59:353-355.
 14. LaFortune S. The use of confirming x-rays to verify tip position for peripherally inserted catheters. *JIN.* 1993;16:246-250.
 15. Loughran SC, Edwards S, McClure S. Peripherally inserted central catheters – guidewire versus non-guidewire use: A comparative study. *JIN.* 1992;15:152-159.
 16. Marx M. The management of the difficult peripherally inserted central venous catheter line removal. *JIN.* 1995;18:246-249.
 17. Maschke SP, Rogove HJ. Cardiac tamponade associated with a multilumen central venous catheter. *Crit Care Med.* 1984;12:611-612.
 18. Peters JL, ed. *A Manual of Central Venous Catheterization and Parenteral Nutrition.* Boston, Mass: John Wright PSG; 1983:140-146, 156.
 19. Ryder MA. Peripherally inserted central venous catheters. *Nursing Clinics of North America.* 1993;28:937-971.
 20. Ryder M. Total parenteral nutrition. In: Boggs RL, Woodridge-Kin M, eds. *AACN Procedural Manual for Critical Care Nursing.* 3rd ed. Philadelphia, PA: WB Saunders Co: 1993:771-812.
 21. Sheep RE, Guiney WB Jr. Fatal cardiac tamponade. *JAMA.* 1982;248:1632-1635.
 22. Sigurdsson J, Riba P, Sigurdsson S. The wandering central venous catheter. *Intensive Care Med.* 1985;11:263-264.
 23. Tocino IM, Watanabe A. Impending catheter perforation of superior vena cava: radiographic recognition. *Am J Roentgenology.* 1986;146:487-490.
 24. Todres ID, Cote CJ. Procedures. In: Cote CJ, Ryan JF, Todres ID, Goudsouzian NG, eds. *A Practice of Anesthesia for Infants and Children.* 2nd ed. Philadelphia, PA: WB Saunders Co; 1993:505-523.
 25. Wall JL, Kierstead VL. Peripherally inserted central catheters resistance to removal: A rare complication. *JIN.* 1995;18:251-254.
- Arrow International, Inc., recommends that the user be acquainted with the reference literature.
- *If you have any questions or would like additional reference information, please contact Arrow International, Inc.
- †Steri-Strip® is a registered trademark of 3M Health Care.

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Arrow International, Inc. rekommenderar att användaren studerar referenslitteraturen.
*Om du har frågor eller önskar ytterligare referensinformation, kontakta Arrow International, Inc.