



ARROW® PICC

ARROW®

REF **CDA-35041-HPK1A**

4 Fr. | **50** cm catheter length | **.018** inch dia. spring-wire guide | Distal Lumen ARROW VPS® Compatible

ErgoPack®


Contents:

- 1: Indwelling TaperFree™ Catheter: 4 Fr. x 50 cm Pressure Injectable Radiopaque Polyurethane with Blue FlexTip®, T-Port Connector, Extension Line Clamps, Contamination Guard and Placement Wire
- 1: GlideThru™ Peel-Away Sheath: 4 Fr. x 2-3/4" (7 cm) Radiopaque over 4 Fr. Dilator
- 1: Spring-Wire Guide, Nitinol, Marked: .018" (0.46 mm) dia. x 17-3/4" (45 cm) (Straight Soft Tip on One End - Straight Stiff Tip on Other) with Arrow Advancer
- 1: Injection Needle: SafetyGlide™ 25 Ga. x 5/8" (1.60 cm) RW
- 1: Safety Introducer Needle: Echogenic 21 Ga. x 2-3/4" (7 cm) TW
- 1: Introducer Needle: Echogenic 21 Ga. x 2-3/4" (7 cm) TW
- 1: Syringe: 10 mL Luer-Lock
- 1: Syringe: 3 mL Luer-Lock
- 1: Dust Cap: Non-Vented
- 1: SecondSite™ Adjustable Hub: Catheter Clamp
- 1: SecondSite™ Adjustable Hub: Fastener
- 1: SharpsAway® Disposal Cup
- 1: SharpsAway II™ Locking Disposal Cup
- 1: Catheter Trimmer
- 1: Towel: 24" x 36" (61 cm x 91 cm)
- 1: Drape: 60" x 76" (152 cm x 193 cm)
- 1: Drape: 36" x 41" (91 cm x 104 cm) with 3" x 5" (8 cm x 13 cm) fenestration
- 1: Towel
- 1: Local Anesthetic Component
- 1: Safety Scalpel: #11
- 1: Flow Rate Information Card with Injection Log
- 1: Sterile Procedure Sign
- 1: Patient ID Card
- 1: Chart Sticker
- 1: Patient Information Booklet
- 1: Checklist/CLIP Sheet
- 2: Paper Tape Measure
- 1: Tourniquet
- 5: Gauze Pad: 2" x 2" (5 cm x 5 cm)
- 5: Gauze Pad: 4" x 4" (10 cm x 10 cm)
- 1: Dressing: Tegaderm® 15.5 cm x 10 cm
- 1: Surgical Apparel: Impervious Gown
- 1: Surgical Apparel: Mask with Eye Shield
- 1: Surgical Apparel: Bouffant Cap
- 1: Surgical Apparel: Mask
- 1: Tape: Steri-Strip®²
- 1: HemoHopper® Fluid Receptacle

¹A trademark of Becton, Dickinson and Company.
²A registered trademark of 3M Company.

All components are CE 0086 unless otherwise noted.

Contraindications: The Pressure Injectable PICC is contraindicated wherever there is presence of device related infections or presence of previous or current thrombosis in the intended insertion vessel or catheter pathway. Clinical assessment of patient must be completed to ensure no contraindications exist.

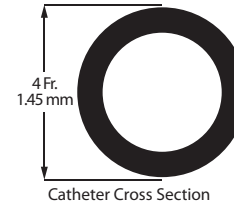

Lumen	Priming Volume* (mL)	Flow Rate† (mL/hr)	Pump Flow Rate‡ (mL/hr)	MAX Pressure Injection Flow Rate** (mL/sec)
Distal (17 Ga.)	0.55	1220	5740	5

* Priming volumes are approximate and are done without accessories.
 † Flow rates were determined using room temperature water, 100 cm head height and represent approximate flow capabilities.
 ‡ Pump flow rates are determined at maximum pump pressure of 10 psig and represent approximate flow capabilities.
 ** Pressure injection flow rates are determined at the injector pressure setting of 300 psi maximum using media of 11.8 centipoise viscosity, with 152 cm pressure tubing.



Not made with natural rubber latex