

Pressure Injectable Two-Lumen CVC Kit

with Blue FlexTip® Catheter



ARROW
INTERNATIONAL

Arrow Select Kits

NA-12702-PCDC

> Sharps Safety Features

Contents:

- 1: Two-Lumen Indwelling Catheter: 7 Fr. x 16 cm Pressure Injectable Radiopaque Polyurethane with Blue FlexTip®, Extension Line Clamps
- 1: Spring-Wire Guide, Marked: .032" (0.81 mm) dia. x 23-5/8" (60 cm) (Straight Soft Tip on One End – "J" Tip on Other) with Arrow Advancer™
- 1: Catheter: 18 Ga. x 2-1/2" (6.35 cm) Radiopaque over 20 Ga. RW Introducer Needle
- 1: Pressure Transduction Probe
- 1: Introducer Needle: Echogenic 18 Ga. x 2-1/2" (6.35 cm) XTW and 5 mL Arrow® Raulerson Spring-Wire Introduction Syringe
- 1: **Injection Needle: SafetyGlide®¹ 25 Ga. x 1" (2.54 cm) and 3 mL Luer-Lock Syringe**
- 1: **Injection Needle: SafetyGlide®¹ 22 Ga. x 1-1/2" (3.81 cm) and 5 mL Luer-Slip Syringe**
- 1: 5 mL 1% Lidocaine HCl Solution
- 1: Tissue Dilator
- 1: **3 mL Applicator 2% CHG and 70% IPA ChloroPrep®² One-Step Solution with Hi-Lite Orange™ Tint**
- 1: Needle Holder
- 1: **Drape: 41" x 55" with 4" fenestration**
- 1: Local Anesthetic Components
- 1: Suture: 3-0 Silk with Curved Needle
- 2: Gauze Pads: 2" x 2"
- 5: Gauze Pads: 4" x 4"
- 1: SecondSite™ Adjustable Hub: Catheter Clamp and Fastener
- 1: **Safety Scalpel: #11**
- 1: **SharpsAway II™ Locking Disposal Cup**
- 1: SharpsAway® Disposal Cup
- 1: Medication Label: 1% Lidocaine
- 1: Flow Rate Information Card with Injection Log
- 1: Tubing: 20" Extension±

Infection Protection Features:

- > ChloroPrep® One-Step skin antiseptic, 2% chlorhexidine gluconate in 70% isopropyl alcohol with Hi-Lite Orange™ Tint
- > Extra-Large Drape size

Safety Features:

- > SafetyGlide® Protected Needles
- > SharpsAway II™ Locking Disposal Cup
- > Safety Scalpel

Rx only.

Warning: Read all package insert warnings, precautions, and instructions prior to use. Failure to do so may result in severe patient injury or death.

All contents sterile unless otherwise indicated.

Single use – do not resterilize. Sterilized by ethylene oxide. Keep away from sunlight. Keep dry. Do not use if package has been previously opened or damaged. Fluid path components are non-pyrogenic. Contains medication.

±**Precaution: Contains Phthalate: DEHP.**


Arrow International provides the enclosed medication label for your convenience. Please ensure that this label is applied to the correct syringe and corresponding medication.

Store between 20 - 25°C (68 - 77°F). Avoid freezing and excessive heat above 40°C (104°F).



Arrow International, Inc.
2400 Bernville Road, Reading, PA 19605 USA
www.arrowintl.com

IC-12702-131B (3/16)

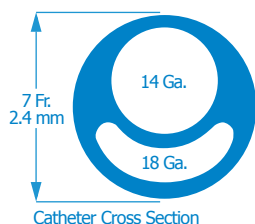

Lumen	Priming Volume* (mL)	Flow Rate† (mL/hr)	MAX Pressure Injection Flow Rate‡ (mL/sec)
Distal (14 Ga.)	0.54	5300	10
Proximal (18 Ga.)	0.37	2100	10

* Priming volumes are approximate and are done without injection caps.
† Based on 100 cm head height.
‡ Pressure injection flow rates are determined using media of 11.8 centipoise viscosity, with 60" administration set tubing.

Pressure Injectable Two-Lumen CVC Kit

NA-12702-PCDC

7 Fr.

2 Lumen

16 cm catheter length

.032 inch diameter spring-wire guide

> Sharps Safety Features

Arrow Select Kits



ARROW
INTERNATIONAL