



# ARROW® PICC powered by ARROW® VPS TipTracker™ Stylet



REF **CDC-35552-TTS**

**5** Fr. | **2** Lumen | **55** cm catheter length | **.018** inch dia. spring-vice guide



**Contents:**

- 1: Two-Lumen Indwelling TaperFree™ Catheter: 5 Fr. x 55 cm Pressure Injectable Radiopaque Polyurethane with Blue FlexTip®, Extension Line Clamps, Contamination Guard, Arrow® VPS TipTracker™ Stylet and T-Port Connector
- 1: GlideThru™ Peel-Away Sheath: 5 Fr. x 2-3/4" (7 cm) Radiopaque over 5 Fr. Dilator
- 1: Placement Wire Assembly
- 1: Spring-Wire Guide, Nitinol, Marked: .018" (0.46 mm) dia. x 17-3/4" (45 cm) (Straight Soft Tip on One End - Straight Stiff Tip on Other) with Arrow Advancer
- 1: Injection Needle: SafetyGlide™ 25 Ga. x 5/8" (1.60 cm) RW
- 1: Safety Introducer Needle: Echogenic 21 Ga. x 2-3/4" (7 cm) TW
- 1: Blunt Fill Needle: 18 Ga. x 1-1/2" (3.81 cm)
- 1: Introducer Needle: Echogenic 21 Ga. x 2-3/4" (7 cm) TW
- 1: Syringe: 10 mL Luer-Lock
- 1: Syringe: 3 mL Luer-Lock
- 1: Skin Protectant Prep Pad (non-sterile solution)
- 4: 10 mL Luer-Lock Pre-Filled Saline Syringe
- 1: 5 mL 1% Lidocaine HCl Solution
- 1: 3 mL Applicator 2% CHG and 70% IPA ChloroPrep®<sup>2</sup> One-Step Solution (non-sterile solution)

- 1: 3 mL Applicator 2% CHG and 70% IPA ChloroPrep®<sup>2</sup> One-Step Solution with Hi-Lite Orange™ Tint (non-sterile solution)
- 2: Dust Cap: Non-Vented
- 1: SecondSite™ Adjustable Hub: Fastener
- 1: SecondSite™ Adjustable Hub: Catheter Clamp
- 1: SharpsAway® Disposal Cup
- 1: SharpsAway II™ Locking Disposal Cup
- 1: Catheter Trimmer
- 1: Towel: 24" x 36" (61 cm x 91 cm)
- 1: Drape: 60" x 76" (152 cm x 193 cm)
- 1: Drape: 36" x 41" (91 cm x 104 cm) with 3" x 5" (8 cm x 13 cm) fenestration
- 1: Towel
- 1: Safety Scalpel: #11
- 1: Flow Rate Information Card with Injection Log
- 1: Sterile Procedure Sign
- 1: Patient ID Card
- 1: Chart Sticker
- 1: Patient Information Booklet
- 1: Medication Label: 1% Lidocaine
- 1: Checklist/CLIP Sheet
- 2: Paper Tape Measure

- 1: Tourniquet
- 1: Foam Electrodes (3 per pack)
- 10: Gauze Pad: 2" x 2" (5 cm x 5 cm)
- 10: Gauze Pad: 4" x 4" (10 cm x 10 cm)
- 1: Dressing: Tegaderm®<sup>3</sup> 15.5 cm x 10 cm
- 1: Surgical Apparel: Impervious Gown
- 1: Dressing: BIOPATCH®<sup>4</sup>
- 1: Dressing: STATLOCK®<sup>5</sup> Catheter Stabilization Device
- 1: Surgical Apparel: Mask with Eye Shield
- 1: Surgical Apparel: Bouffant Cap
- 1: Surgical Apparel: Mask
- 1: Equipment Cover: 13 cm x 244 cm (5" x 96")
- 1: Tape: Steri-Strip®<sup>3</sup>
- 1: HemoHopper® Fluid Receptacle
- 1: T-Piece Cover

<sup>1</sup>A trademark of Becton, Dickinson and Company.  
<sup>2</sup>A registered trademark of CareFusion or one of its subsidiaries.  
<sup>3</sup>A registered trademark of 3M Company.  
<sup>4</sup>A registered trademark of Johnson & Johnson Corporation.  
<sup>5</sup>A registered trademark of C. R. Bard, Inc.

Sterilized using ethylene oxide unless otherwise indicated in the contents list.

**Contraindications: The Pressure Injectable PICC is contraindicated wherever there is presence of device related infections or presence of previous or current thrombosis in the intended insertion vessel or catheter pathway. Clinical assessment of patient must be completed to ensure no contraindications exist.**


Arrow International provides the enclosed medication label(s) for your convenience. Please ensure that these labels are applied to the correct syringe and corresponding medication.



Not made with natural rubber latex



Store between 20-25°C (68-77°F).

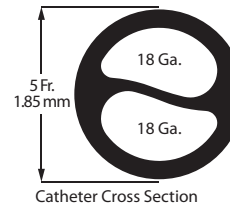

Lumen	Priming Volume* (mL)	Flow Rate† (mL/hr)	Pump Flow Rate‡ (mL/hr)	MAX Pressure Injection Flow Rate** (mL/sec)
Distal (18 Ga.)	0.52	630	3660	5
Proximal (18 Ga.)	0.53	620	3560	5

\* Priming volumes are approximate and are done without accessories.  
 † Flow rates were determined using room temperature water, 100 cm head height and represent approximate flow capabilities.  
 ‡ Pump flow rates are determined at maximum pump pressure of 10 psig and represent approximate flow capabilities.  
 \*\* Pressure injection flow rates are determined at the injector pressure setting of 300 psi maximum using media of 11.8 centipoise viscosity, with 152 cm pressure tubing.