



# Two-Lumen Central Venous Catheterization Kit with Blue FlexTip® Catheter

REF **AK-12702-CDC**

**7** Fr. | **2** Lumen | **16** cm catheter length | **.025** inch dia. spring-wire guide

**Contents:**

- 1: Two-Lumen Indwelling Catheter: 7 Fr. x 16 cm Radiopaque Polyurethane with Blue FlexTip®, Extension Line Clamps
- 1: Spring-Wire Guide, Marked: .025" (0.64 mm) dia. x 17-3/4" (45 cm) (Straight Soft Tip on One End - "J" Tip on Other) with Arrow Advancer
- 1: Catheter: 20 Ga. x 1-3/4" (4.45 cm) Radiopaque over 22 Ga. TW Introducer Needle
- 1: Injection Needle: 25 Ga. x 1" (2.54 cm)
- 1: Introducer Needle: 18 Ga. x 2-1/2" (6.35 cm) RW
- 1: Local Anesthetic Components
- 1: Pressure Transduction Probe
- 1: Introducer Needle: 20 Ga. x 1-1/2" (3.81 cm) XTW and 5 mL Arrow® Raulerson Spring-Wire Introduction Syringe
- 1: Injection Needle: 22 Ga. x 1-1/2" (3.81 cm) and 5 mL Luer-Slip Syringe
- 1: Syringe: 3 mL Luer-Lock

- 1: Tissue Dilator
- 1: 5 mL 1% Lidocaine HCl Solution
- 1: 3 mL Applicator 2% CHG and 70% IPA Chloraprep®<sup>1</sup> One-Step Solution with Hi-Lite Orange™ Tint
- 2: Dust Cap: Non-Vented
- 1: SecondSite™ Adjustable Hub: Fastener
- 1: SecondSite™ Adjustable Hub: Catheter Clamp
- 1: SharpsAway II™ Locking Disposal Cup
- 1: Drape: 41" x 55" with 4" fenestration
- 1: CSR Wrap
- 1: Safety Scalpel: #11
- 2: Gauze Pad: 2" x 2" (5 cm x 5 cm)
- 5: Gauze Pad: 4" x 4" (10 cm x 10 cm)
- 1: Suture: 3-0 Braided Silk with Straight Needle

<sup>1</sup>A registered trademark of CareFusion or one of its subsidiaries.


DuPont™ Tyvek®  
packaging

Rx only. Warning: Read all package insert warnings, precautions and instructions prior to use. Failure to do so may result in severe patient injury or death.

Warning: This product contains a chemical known to the State of California to cause cancer, birth defects or other reproductive harm.

All contents sterile unless otherwise indicated. Fluid path components are non-pyrogenic.

|         |                              |              |                    |                              |                         |          |                                  |                                    |                                      |
|---------|------------------------------|--------------|--------------------|------------------------------|-------------------------|----------|----------------------------------|------------------------------------|--------------------------------------|
|         |                              |              |                    |                              |                         |          |                                  |                                    |                                      |
| Caution | Consult instructions for use | Do not reuse | Do not resterilize | Sterilized by ethylene oxide | Keep away from sunlight | Keep dry | Do not use if package is damaged | Not made with natural rubber latex | Store between 20 - 25°C (68 - 77°F). |

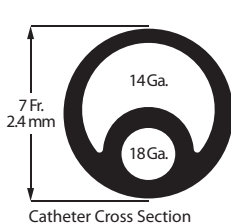

| Lumen             | Priming Volume* (mL) | Gravity Flow Rate† (mL/hr.) |
|-------------------|----------------------|-----------------------------|
| Distal (18 Ga.)   | 0.32                 | 1500                        |
| Proximal (14 Ga.) | 0.47                 | 5000                        |

\* Priming volumes are approximate and are done without accessories.  
† Flow rates are done with deionized water or normal saline, room temperature, 100 cm head height and represent approximate flow capabilities.

Manufacturer:  
Arrow International, Inc.  
Subsidiary of Teleflex Incorporated  
2400 Bernville Road  
Reading, PA 19605 USA

**Teleflex®**  
www.teleflex.com

IK-12702-105F (4/17)



# Two-Lumen Central Venous Catheterization Kit with Blue FlexTip® Catheter

REF **AK-12702-CDC**

**7** Fr. | **2** Lumen | **16** cm catheter length | **.025** inch dia. spring-wire guide