



Arrowg+ard Blue® Multi-Lumen CVC Super Kit

7 Fr. | **3** Lumen | **30** cm catheter length | **.032** inch dia. spring-wire guide

Contents:

- 1: Multi-Lumen Indwelling Catheter: 7 Fr. x 30 cm Radiopaque Polyurethane with Blue FlexTip®, Arrowg+ard® Antimicrobial Surface Treatment¹, Extension Line Clamps
- 1: Spring-Wire Guide, Marked: .032" (0.81 mm) dia. x 26-3/4" (68 cm) (Straight Soft Tip on One End - "J" Tip on Other) with Arrow Advancer
- 1: Catheter: 18 Ga. x 2-1/2" (6.35 cm) Radiopaque over 20 Ga. RW Introducer Needle
- 1: Injection Needle: 25 Ga. x 1" (2.54 cm)
- 1: Blunt Fill Needle: 18 Ga. x 1-1/2" (3.81 cm)
- 1: Pressure Transduction Probe
- 1: Introducer Needle: 18 Ga. x 2-1/2" (6.35 cm) XTW and 5 mL Arrow® Raulerson Spring-Wire Introduction Syringe
- 1: Injection Needle: 22 Ga. x 1-1/2" (3.81 cm) and 5 mL Luer-Slip Syringe
- 1: Syringe: 3 mL Luer-Lock
- 1: Tissue Dilator: 8.5 Fr. (2.8 mm) x 10.2 cm
- 1: 5 mL 1% Lidocaine HCl Solution and Alcohol Prep
- 1: 3 mL Applicator 2% CHG and 70% IPA Chloraprep® One-Step Solution with Hi-Lite Orange™ Tint (non-sterile solution)
- 3: Dust Cap: Non-Vented
- 1: SecondSite™ Adjustable Hub: Fastener
- 1: SecondSite™ Adjustable Hub: Catheter Clamp
- 1: SharpsAway II™ Locking Disposal Cup
- 1: SharpsAway® Disposal Cup
- 1: Drape: 41" x 50" with 4" fenestration, with adhesive
- 1: Safety Scalpel: #11

- 2: Gauze Pad: 2" x 2" (5 cm x 5 cm)
- 5: Gauze Pad: 4" x 4" (10 cm x 10 cm)
- 1: Suture: 3-0 Braided Silk with Straight Needle

¹Licensed under US Patent No. 6,706,024.

²A registered trademark of CareFusion or one of its subsidiaries.

Rx only. **Warning: Read all package insert warnings, cautions, and instructions prior to use. Failure to do so may result in severe patient injury or death.**

Sterilized using ethylene oxide unless otherwise indicated in the contents list.

California Prop. 65

Warning: Cancer and Reproductive Harm www.P65Warnings.ca.gov.

Fluid path components are non-pyrogenic.

Contraindications: The Arrowg+ard Blue antimicrobial catheter is contraindicated for patients with known hypersensitivity to chlorhexidine, silver sulfadiazine and/or sulfa drugs (refer to product instructions for references).

Lumen	Priming Volume* (mL)	Gravity Flow Rate† (mL/hr)
Distal (16 Ga.)	0.55	2400
Medial (18 Ga.)	0.45	1070
Proximal (18 Ga.)	0.50	1160


* Priming volumes are approximate and are done without accessories.
 † Flow rates are done with deionized water or normal saline, room temperature, 100 cm head height and represent approximate flow capabilities.



STERILE EO



Not made with natural rubber latex

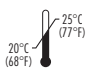

Store between 20-25°C (68-77°F).



Arrowg+ard Blue® Multi-Lumen CVC Super Kit

7 Fr. | **3** Lumen | **30** cm catheter length | **.032** inch dia. spring-wire guide

Manufacturer:
 Arrow International, Inc.
 Subsidiary of Teleflex Incorporated
 2400 Bernville Road
 Reading, PA 19605 USA

Teleflex®
www.teleflex.com

IK-24703-115A (2020-03)