

Multi-Lumen CVC Kit



REF ASK-42703-PI

7 Fr. | 3 Lumen | 16 cm catheter length | .032 inch dia. spring-wire guide

Contents:

- 1: Multi-Lumen Indwelling Catheter: 7 Fr. x 16 cm Pressure Injectable Radiopaque Polyurethane with Blue FlexTip®, Arrowgard Blue PLUS® Antimicrobial Surface Treatment¹, Extension Line Clamps
- 1: Spring-Wire Guide, Marked: .032" (0.81 mm) dia. x 17-3/4" (45 cm) (Straight Soft Tip on One End - "J" Tip on Other) with Arrow Advancer
- 1: Catheter: 18 Ga. x 2-1/2" (6.35 cm) Radiopaque over 20 Ga. RW Introducer Needle
- 1: Pressure Transduction Probe
- 1: Introducer Needle: Echogenic 18 Ga. x 2-1/2" (6.35 cm) XTW and 5 mL Arrow® Raulerson Spring-Wire Introduction Syringe
- 1: Injection Needle: Eclipse™² 22 Ga. x 1-1/2" (3.81 cm) RW and 5 mL Luer-Slip Syringe
- 1: Injection Needle: Eclipse™² 25 Ga. x 1" (2.54 cm) RW and 3 mL Luer-Lock Syringe
- 1: Tissue Dilator: 8.5 Fr. (2.8 mm) x 10.2 cm
- 1: SecondSite™ Adjustable Hub: Fastener
- 1: SecondSite™ Adjustable Hub: Catheter Clamp
- 1: SharpsAway® Disposal Cup
- 1: SharpsAway II™ Locking Disposal Cup
- 1: Drape: 21" x 36" (53 cm x 91 cm) with 3" (7.6 cm) fenestration, with adhesive
- 1: Safety Scalpel: #11
- 1: Flow Rate Information Card with Injection Log
- 2: Gauze Pad: 2" x 2" (5 cm x 5 cm)
- 5: Gauze Pad: 4" x 4" (10 cm x 10 cm)
- 1: Suture: 3-0 Braided Silk with Straight Needle

¹ Licensed under US Patent Nos. 6,706,024; 6,872,195 and 7,329,412.
² A trademark of Becton, Dickinson and Company.




Rx only. **Warning:** Read all package insert warnings, cautions, and instructions prior to use. Failure to do so may result in severe patient injury or death.

California Prop. 65
 Warning: Cancer and Reproductive Harm www.P65Warnings.ca.gov.

Contraindications: The Arrowgard Blue Plus antimicrobial catheter is contraindicated for patients with known hypersensitivity to chlorhexidine, silver sulfadiazine and/or sulfa drugs (refer to product instructions for references).

Sterilized using ethylene oxide unless otherwise indicated in the contents list. Fluid path components are non-pyrogenic.

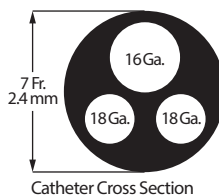


Lumen	Priming Volume* (mL)	Flow Rate† (mL/hr)	MAX Pressure Injection Flow Rate‡ (mL/sec)
Distal (16 Ga.)	0.38	3000	10
Medial (18 Ga.)	0.37	1400	5
Proximal (18 Ga.)	0.39	1600	5

* Priming volumes are approximate and are done without an injection cap.
 † Based on 100 cm head height.
 ‡ Pressure injection flow rates are determined using media of 11.8 centipoise viscosity, with 152 cm administration set tubing.

Manufacturer:
 Arrow International, Inc.
 Subsidiary of Teleflex, Incorporated
 2400 Bernville Road
 Reading, PA 19605 USA



IP-42703-111A (6/19)

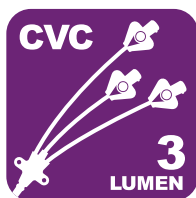

Multi-Lumen CVC Kit



REF ASK-42703-PI

7 Fr. | 3 Lumen | 16 cm catheter length | .032 inch dia. spring-wire guide