



# Peripherally Inserted Central Catheter Set with Blue FlexTip® Catheter Seldinger Access Set with 130 cm Marked Wire

**REF PR-05041-LW**

**4** Fr. | **50** cm catheter length | **.018** inch dia. spring-wire guide

**Contents:**

- 1: Indwelling Catheter: 4 Fr. x 50 cm Radiopaque Polyurethane with Blue FlexTip®, Integral Extension Line, Extension Line Clamp, Integral Suture Wing
- 1: GlideThru™ Peel-Away Sheath: 4 Fr. x 4" (10 cm) Radiopaque over 4 Fr. Dilator
- 1: Spring-Wire Guide: .018" (0.46 mm) dia. x 17-3/4" (45 cm) (Straight Soft Tip on Both Ends)
- 1: Spring-Wire Guide, Marked: .018" (0.46 mm) dia. x 51-3/16" (130 cm)(Straight Soft Tip on One End - Straight Stiff Tip on Other)
- 1: Introducer Needle: Echogenic 21 Ga. x 2-3/4" (7 cm) TW
- 1: Syringe: 5 mL Luer-Slip
- 1: Skin Protectant Prep Pad (non-sterile solution)
- 1: Dust Cap: Non-Vented
- 1: SecondSite™ Adjustable Hub: Catheter Clamp
- 1: SecondSite™ Adjustable Hub: Fastener
- 1: Patient ID Card
- 1: Chart Sticker
- 1: Patient Information Booklet
- 1: Paper Tape Measure
- 1: Dressing: STATLOCK®<sup>1</sup> Catheter Stabilization Device

<sup>1</sup>A registered trademark of C. R. Bard, Inc.

Rx only. **Warning: Read all package insert warnings, cautions and instructions prior to use. Failure to do so may result in severe patient injury or death.**


Sterilized using ethylene oxide unless otherwise indicated in the contents list.

**California Prop. 65**

**Warning: Cancer and Reproductive Harm**  
[www.P65Warnings.ca.gov](http://www.P65Warnings.ca.gov).

Fluid path components are non-pyrogenic.

**Contraindications: The PICC is contraindicated wherever there is presence of device related infections or presence of previous or current thrombosis in the intended insertion vessel or catheter pathway. Clinical assessment of patient must be completed to ensure no contraindications exist.**

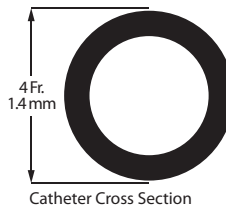

Lumen	Priming Volume* (mL)	Flow Rate† (mL/hr)	Pump Flow Rate‡ (mL/hr)
Distal (18 Ga.)	0.50	850	3700

\* Priming volumes are approximate and are done without accessories.  
 † Flow rates are done with normal saline, room temperature, 100 cm head height and represent approximate flow capabilities.  
 ‡ Pump flow rates are determined at maximum pump pressure of 10 psig.



Arrow International, Inc.  
 Subsidiary of Teleflex Incorporated  
 2400 Bernville Road  
 Reading, PA 19605 USA

**Teleflex®**  
[www.teleflex.com](http://www.teleflex.com)  
 IR-05041-103M (2020-11)



**REF PR-05041-LW**

**4** Fr. | **50** cm catheter length | **.018** inch dia. spring-wire guide

# Peripherally Inserted Central Catheter Set with Blue FlexTip® Catheter Seldinger Access Set with 130 cm Marked Wire