



Pressure Injectable PICC Set with Blue FlexTip® Catheter

and 80 cm Platinum Tipped, Stainless Steel Marked Wire


Distal Lumen Arrow® VPS™ compatible

REF PR-34063-HPSM

Contents:

- ① 1: Multi-Lumen Indwelling Catheter: 6 Fr. x 40 cm Radiopaque **Pressure Injectable** Polyurethane with Blue FlexTip®, Extension Line Clamps and Protection Sleeve
- ② 1: Peel-Away Sheath: 13 Ga. x 4" (10 cm) Radiopaque over 6.5 Fr. Dilator
- ③ 1: Spring-Wire Guide, Marked: .018" (0.46 mm) dia. x 31-1/2" (80 cm) (Straight Soft Tip on One End - Straight Stiff Tip on Other)
- ④ 1: Spring-Wire Guide: .018" (0.46 mm) dia. x 17-3/4" (45 cm) (Straight Soft Tip on Both Ends)
- ⑤ 1: Introducer Needle: Echogenic 21 Ga. x 2-3/4" (7 cm) TW
- ⑥ 1: Syringe: 5 mL Luer-Lock
1: Foil-Laminated Pouch 79% Skin Protectant Prep Pad
3: Dust Caps: Non-Vented
- ⑦ 1: Clamp: Catheter
- ⑧ 1: Fastener: Catheter Clamp
1: Chart Sticker
1: Checklist/CLIP Sheet
1: Patient ID Card
1: Patient Information Booklet
- ⑨ 1: Dressing: STATLOCK®
Catheter Stabilization Device
1: Paper Tape Measure
1: Flow Rate Information Card with Injection Log

¹ A registered trademark of C. R. Bard, Inc.

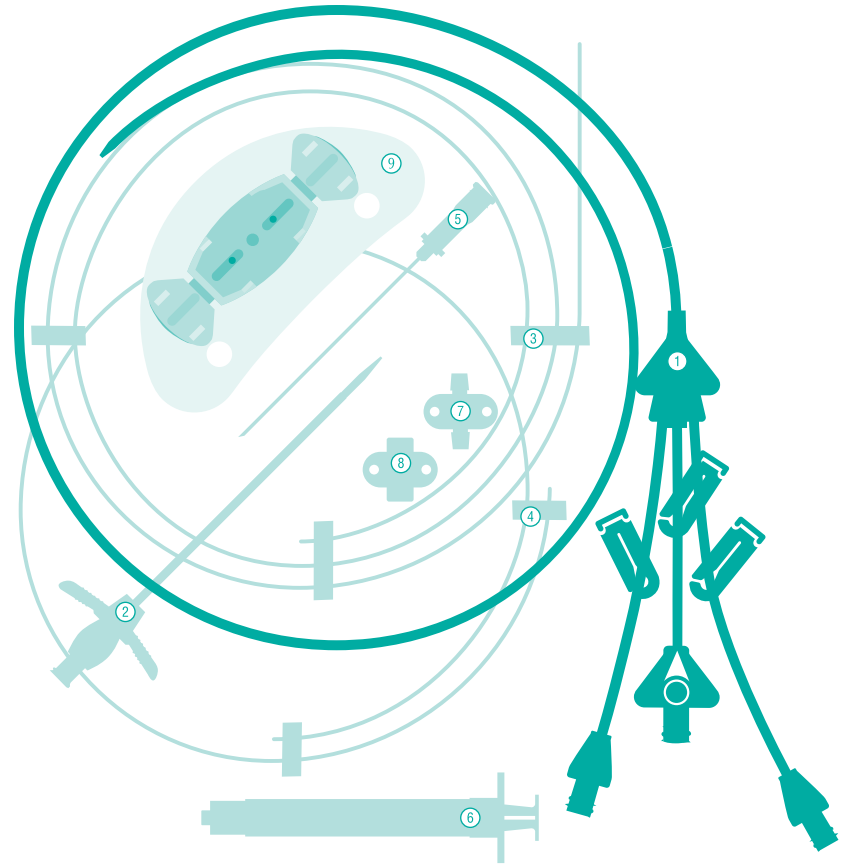

Caution

Do not reuse


Do not resterilize

Sterilized by ethylene oxide

Do not use if package is damaged

Does not contain natural rubber latex

Consult instructions for use

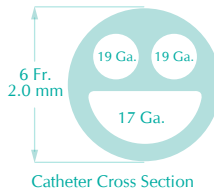

Lumen	Priming Volume* (mL)	Flow Rate† (mL/hr)	Pump Flow Rate‡ (mL/hr)	MAX Pressure Injectable Flow Rate** (mL/sec)
Distal (17 Ga.)	0.60	1650	7526	6
Medial (19 Ga.)	0.40	550	2736	not rated
Proximal (19 Ga.)	0.43	550	2760	not rated

* Priming volumes are approximate and are done without injection caps.
 † Flow rates were determined using room temperature water, 100 cm head height and represent approximate flow capabilities.
 ‡ Pump flow rates are determined at maximum pump pressure of 10 psig and represent approximate flow capabilities.
 ** Pressure injection flow rates are determined at the injector pressure setting of 300 psi maximum using media of 11.8 centipoise viscosity, with 152 cm pressure tubing.

Rx only. Warning: Read all package insert warnings, precautions, and instructions prior to use. Failure to do so may result in severe patient injury or death.

Sterile. Fluid path components are non-pyrogenic. Contains no medication.

ARROW INTERNATIONAL
 2400 Bernville Road
 Reading, PA 19605 USA
www.arrowintl.com

IR-34063-100C (6/16)



Pressure Injectable PICC Set

6 Fr. | 3 Lumen | 40 cm catheter length | .018 inch dia. spring-wire guide

Distal Lumen Arrow® VPS™ compatible

REF PR-34063-HPSM

ARROW INTERNATIONAL

Peel to OPEN