



Pressure Injectable PICC Set with Blue FlexTip® Catheter

and 80 cm Platinum Tipped,
Stainless Steel Marked Wire


Distal Lumen
Arrow® VPS™
compatible

REF PR-35063-HPSM

Contents:

- ① 1: Multi-Lumen Indwelling Catheter: 6 Fr. x 50 cm Radiopaque **Pressure Injectable** Polyurethane with Blue FlexTip®, Extension Line Clamps and Protection Sleeve
- ② 1: Peel-Away Sheath: 13 Ga. x 4" (10 cm) Radiopaque over 6.5 Fr. Dilator
- ③ 1: Spring-Wire Guide, Marked: .018" (0.46 mm) dia. x 31-1/2" (80 cm) (Straight Soft Tip on One End - Straight Stiff Tip on Other)
- ④ 1: Spring-Wire Guide: .018" (0.46 mm) dia. x 17-3/4" (45 cm) (Straight Soft Tip on Both Ends)
- ⑤ 1: Introducer Needle: Echogenic 21 Ga. x 2-3/4" (7 cm) TW
- ⑥ 1: Syringe: 5 mL Luer-Lock
1: Foil-Laminated Pouch 79% Skin Protectant Prep Pad
3: Dust Caps: Non-Vented
- ⑦ 1: Clamp: Catheter
- ⑧ 1: Fastener: Catheter Clamp
1: Chart Sticker
1: Checklist/CLIP Sheet
1: Patient ID Card
1: Patient Information Booklet
- ⑨ 1: Dressing: STATLOCK®¹ Catheter Stabilization Device
1: Paper Tape Measure
1: Flow Rate Information Card with Injection Log

¹ A registered trademark of C. R. Bard, Inc.

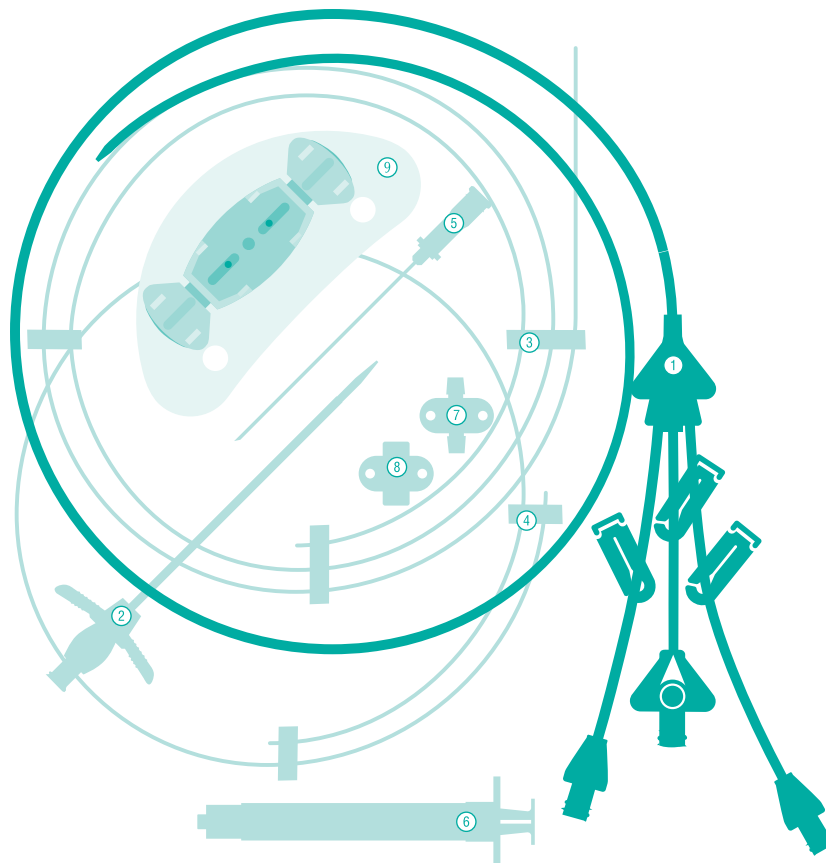

Caution



Do not reuse



Do not resterilize

STERILE EO

Sterilized by ethylene oxide




Do not use if package is damaged



Does not contain natural rubber latex



Consult instructions for use

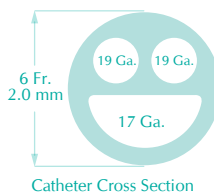

Catheter Cross Section

Lumen	Priming Volume* (mL)	Flow Rate† (mL/hr)	Pump Flow Rate‡ (mL/hr)	MAX Pressure Injectable Flow Rate** (mL/sec)
Distal (17 Ga.)	0.65	1320	6765	6
Medial (19 Ga.)	0.45	450	2597	not rated
Proximal (19 Ga.)	0.46	480	2579	not rated

* Priming volumes are approximate and are done without injection caps.

† Flow rates were determined using room temperature water, 100 cm head height and represent approximate flow capabilities.

‡ Pump flow rates are determined at maximum pump pressure of 10 psig and represent approximate flow capabilities.

**Pressure injection flow rates are determined at the injector pressure setting of 300 psi maximum using media of 11.8 centipoise viscosity, with 152 cm pressure tubing.

Rx only. Warning: Read all package insert warnings, precautions, and instructions prior to use. Failure to do so may result in severe patient injury or death.

Sterile. Fluid path components are non-pyrogenic. Contains no medication.

ARROW
INTERNATIONAL

2400 Bernville Road
Reading, PA 19605 USA
www.arrowintl.com

IR-35063-100C (6/16)



Pressure Injectable PICC Set

6 Fr.

3 Lumen

50 cm catheter length

.018 inch dia. spring-wire guide

Distal Lumen
Arrow® VPS™
compatible

REF PR-35063-HPSM

ARROW
INTERNATIONAL

Peel to OPEN