



ARROW® PICC Set

6 Fr. | 3 Lumen | 55 cm catheter length | .018 inch dia. guide wire

80 cm Hydrophilic-Coated Nitinol Guidewire

Distal Lumen ARROW® VPS® Compatible


REF PR-35563-HPHNM

Contents:

- ① 1: Multi-Lumen Indwelling TaperFree™ Catheter: 6 Fr. x 55 cm **Pressure Injectable** Radiopaque Polyurethane with Blue FlexTip®, Extension Line Clamps, and Contamination Guard
- ② 1: GlideThru™ Peel-Away Sheath: 6 Fr. x 4" (10 cm) Radiopaque over 6 Fr. Dilator
- ③ 1: **Guidewire, Marked: .018" (0.46 mm) dia. x 80 cm Hydrophilic-Coated Nitinol with Soft Tip Tungsten Coil** (CE 0120)
- ④ 1: Guidewire, Marked: .018" (0.46 mm) dia. x 17-3/4" (45 cm) Nitinol (Straight Soft Tip on One End – Straight Stiff Tip on Other)
- ⑤ 1: Introducer Needle: Echogenic 21 Ga. x 2-3/4" (7 cm) TW
- ⑥ 1: Syringe: 10 mL Luer-Lock
- ⑥ 1: Skin Protectant Prep Pad
- ③ 3: Dust Caps: Non-Vented
- ⑦ 1: SecondSite™ Adjustable Hub: Catheter Clamp and Fastener
- ⑧ 1: **Catheter Trimmer**
- ⑨ 1: **Safety Scalpel: #11**
- ① 1: Chart Sticker
- ① 1: Checklist/CLIP Sheet
- ① 1: Flow Rate Information Card with Injection Log
- ① 1: Patient ID Card
- ① 1: Patient Information Booklet
- ⑩ 1: **Dressing: STATLOCK®¹ Catheter Stabilization Device**
- ① 1: Paper Tape Measure


¹ A registered trademark of C. R. Bard, Inc.

All components are CE 0086 unless otherwise noted.

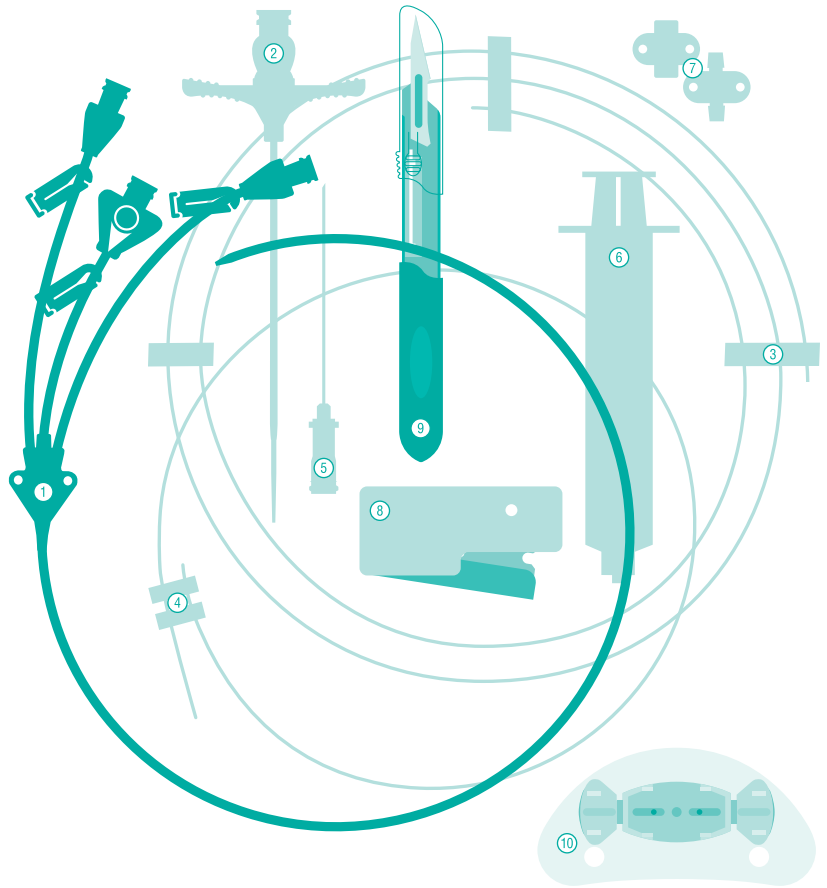

Rx only. Warning: Read all package insert warnings, precautions, and instructions prior to use. Failure to do so may result in severe patient injury or death.

Sterile. Fluid path components are non-pyrogenic. Contains no medication.

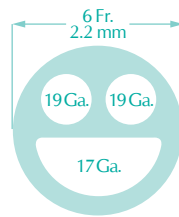

Catheter Cross Section

Lumen	Priming Volume* (mL)	Flow Rate† (mL/hr)	Pump Flow Rate‡ (mL/hr)	MAX Pressure Injection Flow Rate** (mL/sec)
Distal (17 Ga.)	0.71	1070	5976	6
Medial (19 Ga.)	0.46	360	2246	not rated
Proximal (19 Ga.)	0.49	390	2467	not rated

* Priming volumes are approximate and are done without injection caps.

† Flow rates were determined using room temperature water, 100 cm head height and represent approximate flow capabilities.

‡ Pump flow rates are determined at maximum pump pressure of 10 psig and represent approximate flow capabilities.

** Pressure injection flow rates are determined at the injector pressure setting of 300 psi maximum using media of 11.8 centipoise viscosity, with 152 cm pressure tubing.

ARROW
INTERNATIONAL
2400 Bernville Road
Reading, PA 19605 USA
www.arrowintl.com

EC REP Teleflex Medical
IDA Business and
Technology Park
Dublin Road, Athlone
Co. Westmeath, Ireland